

# **CORIO CD-B19 Heating Circulator with Open Bath**


# The new generation of laboratory circulators for routine work and standard tasks

The Open Heating Bath Circulators of the new CORIO™ series feature high-quality bath tanks made of stainless steel and pump connections.

# Your advantages

- · Models for internal and external applications
- Bright, white, easy to read display
- Very quiet
- Easy pump change-over between internal and external circulation
- USB connection
- · High-quality bath tanks made of stainless steel
- Integrated drain screw (except B5)



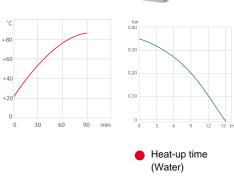

# Pump capacity (Water)

# **Technical Data**

Technical Data	
Order No.	9012419
Model series	CORIO
Category	Heating Circulators with Open Bath
Working temperature range (°C)	+20 +150
Temperature stability (°C)	±0.03
Setting / display resolution	0.01 – 0.1 °C
Temperature Display	LED
Heating capacity (kW)	2
Pump capacity flow rate (I/min)	15
Pump capacity flow pressure (bar)	0.35
Bath opening / bath depth (W x L / D cm)	30 x 35 / 15
Pump connections	M16 x 1
Filling volume liters	14 19
Ambient temperature	540 °C
Dimensions W x L x H (cm)	38 x 58 x 42
Weight (kg)	11.5
Included with each unit	2 barbed fittings suitable for tubing 8 and 12 mm inner dia. (pump connections M16x1 male)
Cooling coil	optional
Bath tank	Stainless steel
Bath cover	optional
Available voltage versions	230 V / 50 Hz 208 V - 230 V / 60 Hz 115 V / 60 Hz 100 V / 50-60 Hz





#### Characteristics



# 100% Checked.

100% testing. 100% quality. Each CORIO™ undergoes thorough quality testing before leaving the factory.



# Modern. Reliable.

High-grade components in every CORIO platinum sensors, proven motor technology, CAN-Bus communication and much more.



#### ATC.

Absolute Temperature Calibration, 1-point calibration (CD).



#### Quick start.

Individual JULABO consultation and comprehensive manuals at your disposal.



#### Brilliant.

Bright display with strong illumination can be viewed easily from a long distance.



#### Safety.

CORIO CD complies with Class III (FL) according to DIN 12876-1 and switches off automatically in case of high temperature or low liquid level alarm.



# Connection. Easy.

Inclined pump connections (M16x1) facilitate the connection of applications. Each unit includes 2 barbed fittings of 8/12 mm diameter each.



#### Satisfied customers.

11 subsidiaries and more than 100 partners worldwide guarantee fast and qualified JULABO support.



# Everything at the front.

All operating controls and safety functions are accessed easily and comfortably from the front.



# Services 24/7.

Around the clock availability. You can find suitable accessories, data sheets, manuals, case studies, and more at www.julabo.com.



# Green technology.

Development consistently applied environmentally friendly materials and technologies.



# Switch on. And off you go.

Intelligent operating concept. Ready for operation with just a few quick and easy steps.



#### Internal. External.

The pump is controlled via a lever located directly below the display. Easily change between internal and external circulation.



#### Timer. Integrated.

CORIO circulators include an integrated timer function. When the set time has elapsed, a signal sounds and the device switches off. Setting range: 0 ... 999 minutes




# JULABO. Quality.

Highest standards of quality for a long product



#### Unique.

Specially designed pyramids ensure that the condensed liquid drips directly back into the bath.

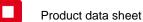





# Locked in.

The lockable power plug guarantees safe connection. More process safety.



# USB.

Remote control made easy. CORIO CD circulators feature a USB connection.